

# Kwawu Arm 3.0

## Socket Version

## Instruction Guide



Jacquín Buchanan

January 2022

This work is licensed under the Creative Commons Attribution 4.0 International License. To view a copy of this license, visit <https://creativecommons.org/licenses/by/4.0/> or send a letter to Creative Commons, PO Box 1866, Mountain View, CA 94042, USA.

In simple language you can:

- Share — copy and redistribute the material in any medium or format.
- Adapt — remix, transform, and build upon the material for any purpose, even commercially.

Under the following terms:

- Attribution — You must give appropriate credit, provide a link to the license, and indicate if changes were made. You may do so in any reasonable manner, but not in any way that suggests the licensor endorses you or your use.
- No additional restrictions — You may not apply legal terms or technological measures that legally restrict others from doing anything the license permits.

This instruction guide shows the steps and processes in making a Kwawu 3.0 Socket Version transradial prosthetic arm and hand. This socket prosthetic arm is best suited for a recipient with a residual forearm that is half or shorter in length compared to the recipient's intact forearm. See the Wrap version of the Kwawu 3.0 for an arm better suited for recipients with a residual arm longer than half the intact arm.

**Improvements of the Kwawu 3.0 over the Kwawu 2.0 device include:**

- Easier stringing the fingers to hand mechanism
- Better response of elbow bend to finger squeeze
- A locking but adjustable wrist
- Improved latch mechanism to fix hand fingers in a grasped position
- The Cuff has a plastic and leather option
- Thermoforming is not required
- A redesigned tensioner
- Easier assembly
- Compatible with wider range of TPU ShoreA hardnesses

## Qualifying Expectations

This arm was designed by hobbyists in hopes that it would fill a need for someone who has had a congenital or traumatic amputation. We do not claim to be professional prosthetists and this arm is designed as an assistive device not a custom prosthetic. This was designed to be used as an assist in everyday functional activities. These are step by step instructions for fabrication of the arm. Our team is available to help problem solve and answer questions should you have difficulties. If you are in need of this assistive device but are unable to make it or have it made we are happy to assist with that as well. As hobbyists our goal is to help make an assistive device which makes day to day activities a little easier.

To complete the instructions in this document and fully prepare an arm for usage can take 30 to 40 hours of labor. Assuming most people reading this document are doing this project in their spare time, expect it to take several weeks. Likely the first one you make may not fit or work exactly right, when you try it on the recipient, so you will have to reprint a part and partially reassemble the arm. It is important to set this expectation for yourself and the person you are making the arm for. Explain to them and yourself that the whole process can take months.

**This Process of making a Kwawu 3.0 Socket Arm is described in the document and is divided into 3 sections, namely:**

**Section 1: Preparing Files for Printing**

**Section 2: 3D Printing of the Necessary Parts**

## Section 3: Assembly of the Kwawu Arm 3.0 Socket Version

# Section 1: Preparing Files for Printing

## Overview


1. Download the Kwawu 3.0 Socket version from [www.thingiverse.com](http://www.thingiverse.com). DO NOT print the stl files directly—they will not result in a functional device—files must be translated into usable stl files using the customizer in OpenSCAD.
2. Download the latest version of OpenSCAD from <http://openscad.org/downloads.html>. Then open the file “Kwawu 3.0- Socket.scad”.
3. Obtain required measurements from the recipient — Record values for the recipient in Table 1.
4. Fill in the OpenSCAD parameters with the specific values of the recipient’s measurements.
5. Render all appropriate parts in the display and save the result as a file.

## Download Kwawu 3.0 files


When you download the files from thingiverse there is one .scad file and a bunch of .STL files. **DO NOT** print the STL files directly. The Open SCAD software will load the STLs and “modify “ them to fit the recipient.

## Download and Install OpenSCAD

1. Download and install the latest version of OpenSCAD. <http://openscad.org/downloads.html>
2. Once you have installed OpenSCAD, click on “New” then in the top menu, select File>Open. Browse and click on the file “Kwawu 3.0 Socket.scad”
3. Uncheck ‘Hide Customizer” in the View menu (as shown below) to show the Customizer UI at the right.

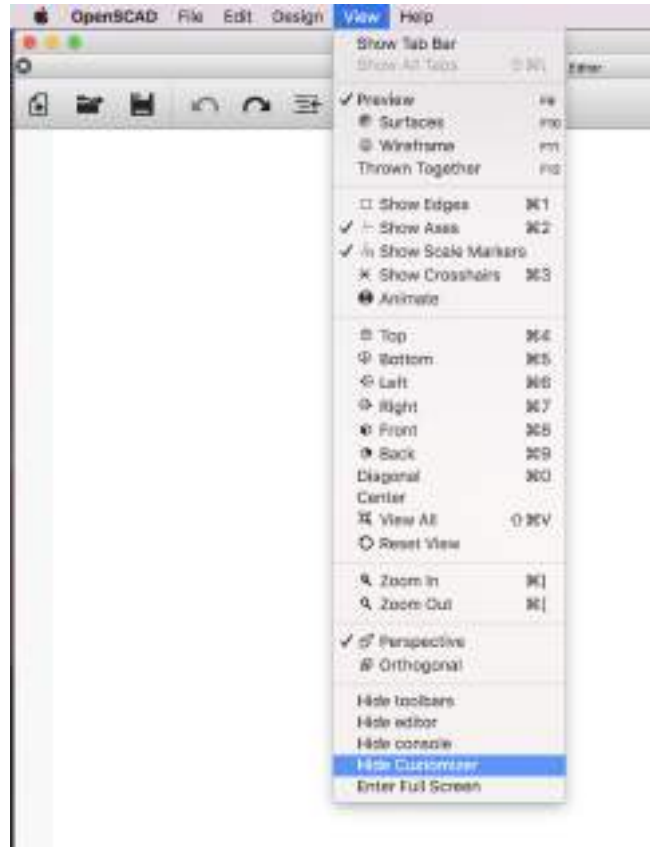


4. In the Customizer menu on the Right Hand side, click the little black arrow to display the parameters. Make sure “Automatic Preview” and “Show Details” are both checked.

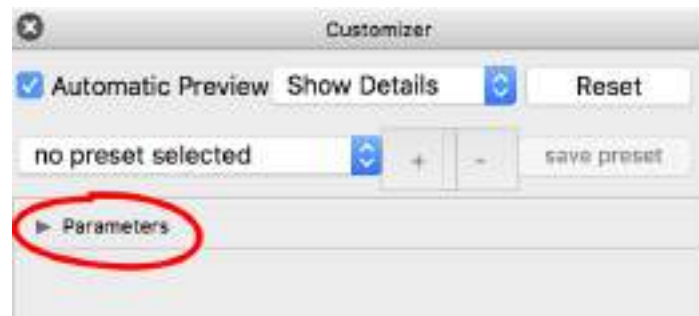

5. Hit the Preview button (F5) to get a quick preview (should not be necessary if ‘Automatic preview’ is enabled).




6. You should see a preview of the Cuff1 in the main window. It looks something like this.

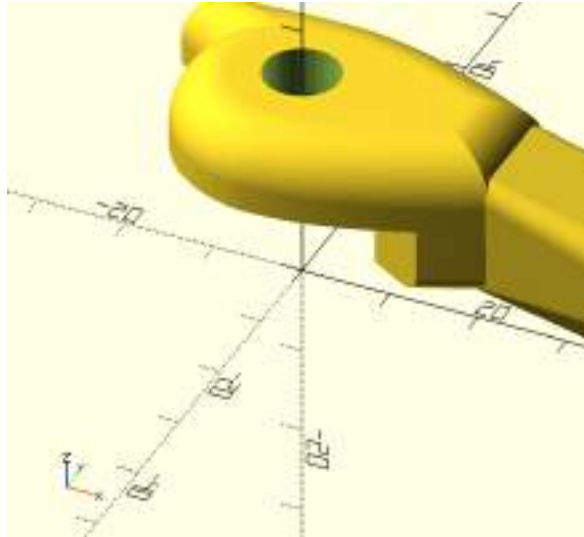

7. If the object is too large or small press the “View all” button




The view will change to something like this.

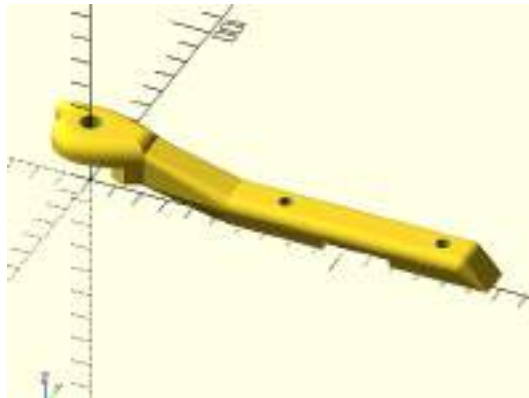

Now you are ready to set all the sizing parameters. Once you have set all the parameters you will need to render each part and save to file.

## Measuring Your Recipient

There are several measurements you will need from the recipient. The circumference measurements you will want to make on the residual arm. The Forearm Length and Hand Width you will make on the intact arm. All measurements are in millimeters.

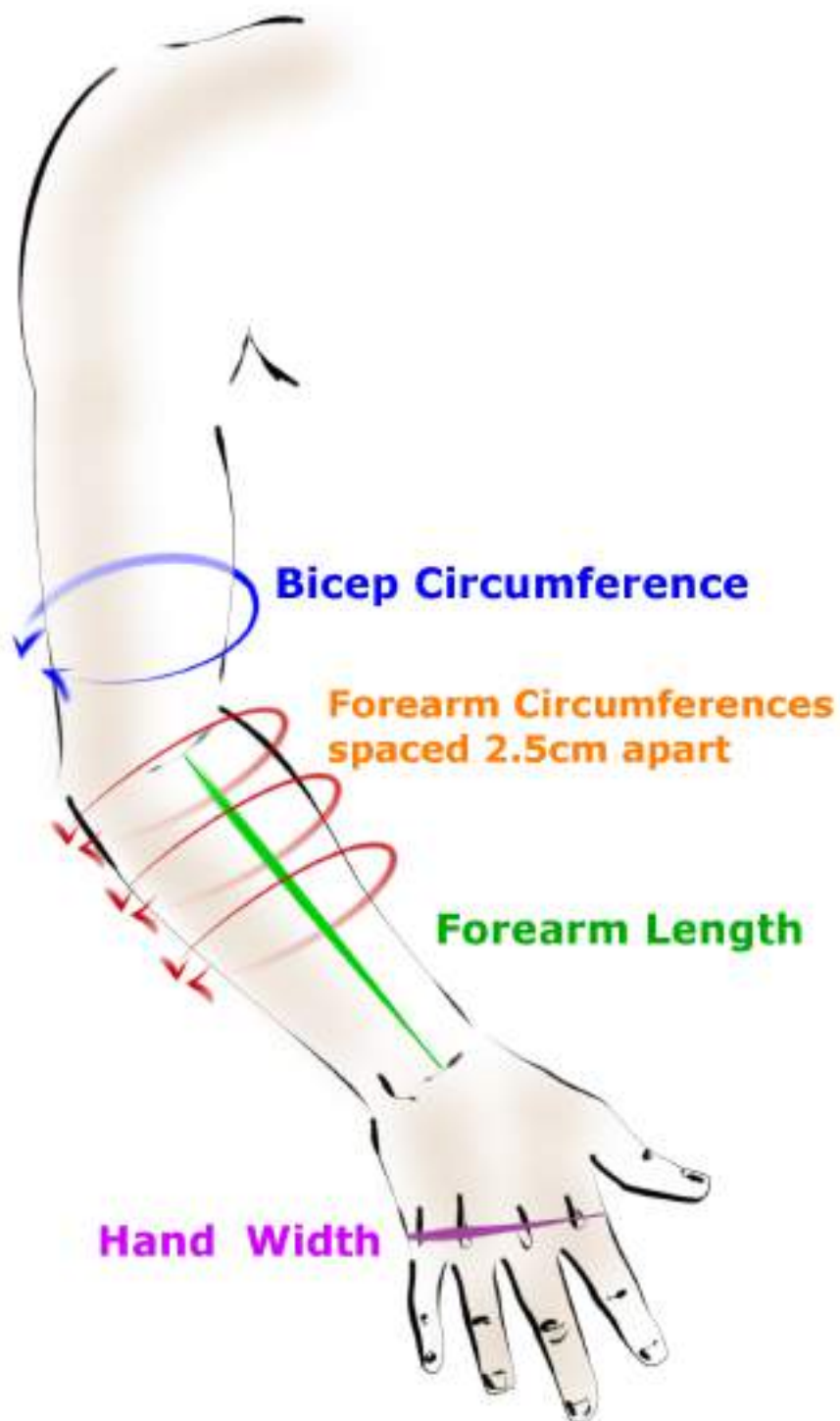

Measurements Diagram

**Bicep Circumference:** Measured around the upper arm at its thickest point. Measured on the residual arm being fitted with the prosthesis.

**Forearm Circumferences:** This is a series of measurements meant to capture the shape of the forearm socket. Starting at the recipient's elbow crease, measure the arm's circumference and every subsequent 25mm, until you reach the end. Measured on the residual arm being fitted with the prosthesis.

**Forearm Length :** Measured from the elbow crease to the wrist crease on the recipient's intact arm. We recommend you use a measurement at least 5% to 10% shorter than the intact arm to improve maneuverability and better mechanical advantage of the prosthetic arm. If the recipient does not mind, you can make the arm length only slightly longer than the residual arm, (leaving enough space for the wrist mechanism)--this will have a less natural look but may be more functional. It is a preference to discuss with the recipient.

**Hand Width:** Measured on the intact arm. The measurement is across the knuckles in millimeters.

## Fill in the OpenSCAD parameters

Once you have measurements from the recipient you will enter various parameters into the Customizer Parameters side bar.

The below example is a screen capture of the OpenSCAD customizer for the socket version. The measurements from the recipient need to be inputted in the appropriate parameters as well as other choices that need to be decided upon consulting the recipient. All these inputs will be discussed separately.

Customizer

Automatic Preview

Show Details

Reset

design default values

+

-

save preset

Parameters

Part

Choose Part

Cuff1

LeftRight

Left or Right Arm

Right

LeatherOrPlastic

Wraps are made from leather or plastic

Leather

HandWidth

Across all four knuckles (mm)

96

ArmLength

Elbow crease to wrist attachment (mm)

271

ArmSplitLength

At what length to split between Upper and Lower arm pieces. Measured from elbow (mm)

135

ForearmCircumferences

Circumferences(mm) of Forearm start at elbow creases separate by 25mm

[272, 268, 265, 261, 255, 0, 0, 0, 0, 0]

BicepCircumference

- Circumference of Bicep (mm)

294

PaddingThickness

Padding Thickness -inside forearm and cuff (mm)

2

PalmBoltDiameter

ISO metric bolt holding palm together (mm)

6

ElbowBoltDiameter

ISO metric bolt holding cuff and arm together (mm)

8

RivetShaftDiameter

Diameter of leather rivet shaft (mm)

5

StrapWidth

Velcro strap width (mm)

25

PencilHolder

Include Pencil Holder

No

UpperVentingHoles

Add Air Venting Holes In Upper Forearm

None

LowerVentingHoles

Add Air Venting Holes In Lower Forearm

None



In OpenSCAD, It may be easier to turn OFF **Automatic Preview** while entering parameters. Many of these parts are slow to render so having the preview run while entering numbers can make the user interface cumbersome.




Detailed information on choices is necessary to render parts to make printable stl files. All 3D parts to be printed must be rendered using the selected choices matching the recipient's measurements and/or preferences. A description of these choices are listed below in order. Table 1 shows the acceptable ranges of these choices and allows you to fill in specific values for your recipient

### OpenSCAD Parameters

Parameter Name	Description
Part	<p>This is a list of all the parts necessary to assemble a Kwawu3 Socket Arm. Choose each of the 33 parts, one by one, then render and save each one.</p> <p>Before rendering any parts, be sure all the other parameters are set to match the recipient. Input or use the sliders to set the values. Uncheck "Automatic Preview" if you have issues with inputting values. After rendering, go to File&gt;Export&gt;As STL . Name the file to correspond to the part rendered and Save. See "Leather or Plastic" for an exception in the export file choice.</p>
Left or Right Arm	The side of the body for the prosthetic arm required by the recipient.
LeatherOrPlastic	<p>The leather option is for the cuff (bicep) area of the device for the socket version. Leather may be a more comfortable option over the plastic version. If the Leather option is chosen, the part "Cuff Leather Template" must be rendered and the result exported: go to File&gt;Export&gt;Export as SVG. Name the file compatible with the template such as "CuffLeatherTemplate " and save. Note that Cuff3 will not render if the Leather Option is chosen.</p>
Handwidth	Measure the width of the recipient's intact hand at the knuckles.
ArmLength	<p>Arm length is typically chosen based on the recipient's "good" arm—or it may be purposely made shorter for easier maneuverability and better mechanical advantage. The socket version is primarily meant for recipients with a residual arm less than half as long as the recipient's intact arm.</p>
ArmSplitLength	<p>If the arm length chosen is too long for your printer, the forearm can be printed in two parts—the UpperArm and LowerArm. If two parts are required, make sure you choose the split point longer than the recipient's residual, and not too close to the wrist (i.e. Do not let the arm split interfere with the wrist bolt). If you choose to make the forearm in one print, set the ArmSplitLength longer than the ArmLength—there is then no need to render LowerArm.</p>

Forearm Circumferences	Insert in the parentheses a series of forearm circumference measurements (comma separated) starting at the recipient's elbow crease and every subsequent 25mm. Once past the recipient's residual length, leave 0's as the remaining numbers. Generally the circumferences should always be decreasing in value, though this is not required. If a circumference step increases, be sure the recipient can still fit their arm in the socket past the tight spot.
Bicep Circumference	Measure the circumference of the bicep where the cuff will wrap around.
Padding Thickness	<p>Padding is inserted in the upper arm area where the recipient's forearm residual will contact the prosthetic. The more comfortable option is silicone sheeting —e.g. <a href="https://amputeestore.com/products/silipos-pressure-relief-gel-padding?variant=4114763009">https://amputeestore.com/products/silipos-pressure-relief-gel-padding?variant=4114763009</a> — but the more expensive option.</p> <p><u>Closed cell neoprene is a less expensive but still a good option. Here are links to a couple of these:</u></p> <p><a href="https://www.amazon.com/Insulation-Neoprene-Adhesive-Multi-Function-Soundproof/dp/B08R9LF7BX">https://www.amazon.com/Insulation-Neoprene-Adhesive-Multi-Function-Soundproof/dp/B08R9LF7BX</a></p> <p><a href="https://www.amazon.com/NATGAI-Neoprene-Adhesive-Multiple-Dimensions/dp/B07VY7SFRP">https://www.amazon.com/NATGAI-Neoprene-Adhesive-Multiple-Dimensions/dp/B07VY7SFRP</a></p> <p><u>2mm is the typical thickness for padding.</u></p>
PalmBolt Diameter	Bolts for the palm can be printed, however, we recommend small bolt sizes (diameter 4mm or less) should be purchased as metal or nylon bolts to reduce the chance of bolt breakage. Also, we recommend the sizes chosen to be in proportion to the arm size, i.e. use small bolts for a small arm and a larger one for larger arms.
ElbowBolt Diameter	Bolts for the elbow can also be printed, however, we recommend small bolt sizes (diameter 4mm or less) should be purchased as metal or nylon bolts to reduce the chance of bolt breakage. Also we recommend the sizes chosen to be in proportion to the arm size, i.e. use small bolts for a small arm and larger one for larger arms.
RivetShaft Diameter	Rivet shaft diameter only applies to the leather version and to how large the holes should be in Cuff 1 and Cuff 2 parts. Rivet shafts come in different sizes. An alternative to rivets are Chicago screw/post sets.
StrapWidth	Velcro Straps are required on the Cuffs for both the leather and plastic versions. Strap widths of 25 mm are typically used for larger sized arms and 10mm for small sizes.

PencilHolder	<p>A Pencil holder on the palm is an option and should be made with TPU.</p> 
UpperVentHoles	<p>Vent holes in Arm1 (UpperVentHoles) are an option for the Kwawu 3 socket version. Different designs or no vent holes are shown so you can choose the pattern (or no venting) as desired. Vent holes are recommended in the area below the arm residual.</p>
LowerVentHoles	<p>Vent holes in the Arm2 (LowerVentHoles) are also an option for the Kwawu 3 socket version. See description in UpperVent Holes.</p>

**Table 1** shows a list of choices that must be made in the OpenSCAD customizer along with ranges and units of each. The column on the right allows the maker to record the choices used for a specific recipient. Many of the choices require measurement on the recipients.

	Table 1: OpenSCAD Options and Ranges		
	Parameter Range	units	Recipient Values
<b>Part</b>	See separate list		n/a
<b>Left or Right</b>	Left or right		
<b>Leather or Plastic</b>	Leather or Plastic		
<b>Handwidth</b>	65 - 186	mm	
<b>ArmLength</b>	141-564	mm	
<b>ArmSplitLength</b>	10-564	mm	
<b>ForearmCircumference</b>	see prior description	mm	
<b>BicepCircumference</b>	160-600	mm	
<b>PaddingThickness</b>	0-10	mm	

<b>PalmBoltDiameter</b>	4, 6, or 8	mm	
<b>ElbowBoltDiameter</b>	4, 6, 8, 10, 12 or 14	mm	
<b>RivetShaftDiameter</b>	2.4, 2.8, 3.1, 4, 5, or 6	mm	
<b>StrapWidth</b>	10-60	mm	
<b>PencilHolder</b>	yes or no		
<b>UpperVentHoles</b>	0, option 1, 2, or 3		
<b>LowerVentHoles</b>	0, option 1, 2, or 3		

## Render all appropriate parts

Select each part in the Parts dropdown. Then Press Render.



NOTE: The render button can take a long time for some parts (**20 minutes or more**).

Save STL (exception—leather template). Export and save the stl files under a name that identifies each part. For the CuffLeatherTemplate choose FILE->EXPORT->EXPORT As SVG



## Section 2: 3D Printing of Necessary Parts



Table 2 shows all parts available for printing along with quantity required, material typically used, and if supports (bed contact) are generally used. PETG is the material recommended for most of the parts because of its strength and solvent resistance—other materials may be tried but we cannot predict a good outcome. Use a suitable TPU for hinges and the pencil holder.

### Printing Recommendations

Typical print temperature and print speed recommended for the material (PETG or TPU).

Infil should be 35% and of honeycomb geometry.


If you are using a printer with a 0.5mm nozzle, you may find a linewidth of 0.4mm works better for some parts (e.g. latch).

We recommend the number of top and bottom layers to be 10 layers. The reason is that the 2mm flexible plastic grid in the case of the plastic cuff has been engineered for flexibility and strength as a solid device. If 10 layers for both top and bottom layers are used, we can be certain the grid will be solid but changes won't be necessary for other parts. Make sure the gaps in the grid are open—if not you may need to raise the z height so the bottom layer/layers will not be spread out which could cause gaps to be fused closed. If gaps are closed, the grid may be inflexible.


**Note: All the parts are rendered from OpenSCAD in the orientation we think is best for printing. Though, you may need to adjust the orientation based on your particular printer.**

A brim is recommended for the finger ends and phalanxes to enhance part stability during printing.

Printing the finger pieces can also benefit from supports/scaffolding. Generally it is better to use a slicing program that will let you manually place the scaffolding for these parts. If you are having trouble with the fingers coming loose from the base plate, try learning them back to a 45 degree angle and adding more scaffolding. This may add to the time to clean the parts, but can be more reliable for some printers.

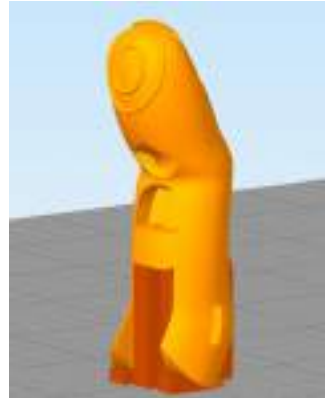

**Note: If scaffolding is used, it is better to be manually added than be auto-generated by the printing software. Auto-generated scaffolding can block holes and make it very difficult to string the hand. Make sure to review the scaffolding so it DOES NOT get created inside a channel or hole where it can be very difficult to remove.**

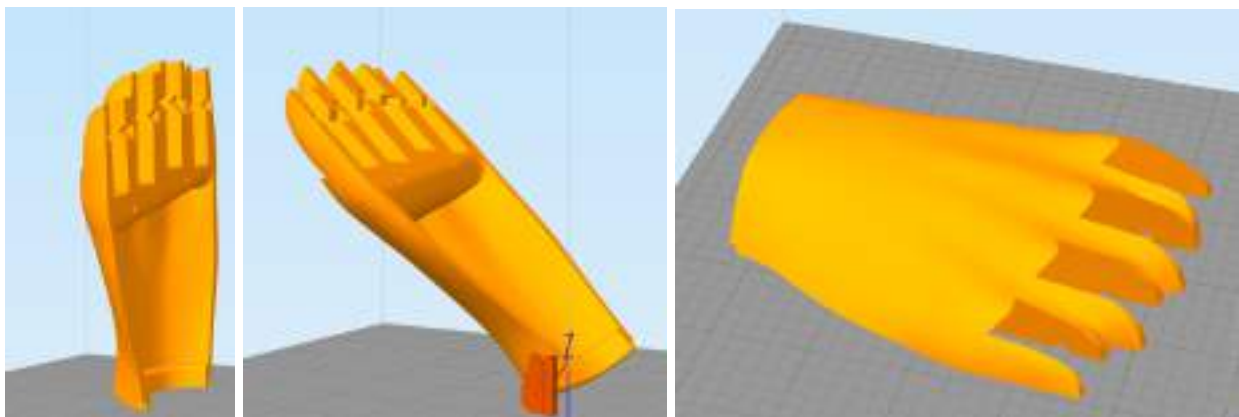

Three orientations to print the Palm Top. The 45 degree angle was suggested by Krzysztof Grandys and is my preferred orientation for best look and strength.

If leather was chosen for the cuff, open the SVG file (saved in the OpenSCAD procedure) in your browser. We have found the **Chrome Browser** opens with an accurate size representation so that you can simply print out the template on paper. We have also used Safari to open the SVG file but found it necessary to size the printing to match the holes in the cuffs. Once an accurate template is printed, lay it out on the leather and cut the outline. Plan to use the smoothest side of the leather in contact with the recipient's skin for greater comfort. A leather punch is recommended to make accurate holes in the leather where it will attach to the Cuff1 and Cuff2 plastic. We also recommend ~4.5 oz leather (~1.8mm thick). It is the more expensive option but we believe it to be more comfortable than the plastic. A website of one supplier of cowhide leather is [www.montanaleather.com](http://www.montanaleather.com)

One choice that is not required in the OpenSCAD options is arm color. The recipient should be consulted about their choice in color for their arm so the printer can purchase a filament color as close as possible to the recipient's choice. If a suitable color is not found, it is possible to use spray paint to better match the recipient's choice. The arm may also be decorated according to the recipient's wishes.

Total printing time for all parts may be 30 to 40 hours. After each print, remove any support structures and sand or file any rough edges or corners. Additional sanding or filing may be necessary when fitting parts together during assembly. Printing several small parts in the same print job will cut down on total print time.

**TPU is required for hinges and the pencil holder. Hinges have been design engineered to use a TPU of Shore A hardness range of 85A to 95A. ShoreA hardness 85A is preferred but more difficult to print than 95A. The brand NinjaFlex is the original material from a company called NinjaTek. They have since come out with several other flexible filaments all with different properties, so be sure to get the one with Shore A Hardness of 85A. Likewise, FilaFlex is a company whose original product was Shore Hardness 82A. They have since come out with FilaFlex Medium which is Shore A Hardness 95A. Shore A Hardness 95A will work for Kwawu 3.0 hinges but since it is stiffer than ShoreA Hardness of lesser values, it would be less responsive.**

**For hinges, make sure the printed parts fit into their appropriate slots snugly, i.e. will not fall out and will fit into the slots without trimming. Since printers can vary, it may be necessary to scale the print in your slicer to obtain a proper fit.**



Table 2: Printing the Parts

Parts	Quantity	Material	Supports	Parts	Quantity	Material	Supports
Palm	1	PETG	yes	ElbowBolt2	0 or 1	PETG/ metal/ nylon	no
PalmTop	1	PETG	optional	Index FingerEnd	1	PETG	optional
WristButton	1	PETG	no	Index Finger Phalanx	1	PETG	optional
WristCover	1	PETG	no	Middle FingerEnd	1	PETG	optional
Cuff1	1	PETG	yes	Middle Finger Phalanx	1	PETG	optional
Cuff2	1	PETG	yes	Pinky FingerEnd	1	PETG	optional
Cuff3	0 or 1	PETG	no	Pinky Finger Phalanx	1	PETG	optional
CuffLeatherTemplat	0 or 1		n/a	RingFingerEnd	1	PETG	optional
UpperArm	1	PETG	optional	RingFingerPhalanx	1	PETG	optional
LowerArm	1	PETG	no	ThumbEnd	1	PETG	optional
Ratchet	1	PETG	yes	Thumb Phalanx	1	PETG	optional
Latch	1	PETG	no	Tensioner	1	PETG	no
Latch Cover	1	PETG	no	Whipple Tree Primary	1	PETG	no
WristBolt	1	PETG	no	Whipple Tree Secondary	1	PETG	no
WristBolt PinCover	1	PETG	no	Pencil Holder	0 or 1	TPU	no
PalmBolt	0 or 2	PETG/ metal/ nylon	no	Hinges	1 set	TPU	no
ElbowBolt1	or	PETG / metal/ nylon	no				



## Section 3: Kwawu 3 Assembly-Socket Version

### How does the Arm Work?

Familiarize yourself with the function of the arm prior to building. The fingers are supposed to close with the thumb first and the pinky last. When the fingers close in order, you get the precision of the thumb and index finger first, while the remaining fingers are still out of the way. This enables you to pinch or to pick up finer objects. If you continue to close the grip the other fingers engage as would happen around a bottle or baseball.

The wrist button allows the wrist to freely rotate when pressed and “lock” in place when released. Different wrist positions are useful for different applications.

The latch allows the grip to stay closed tight with the elbow moving freely. With the latch open, the grip opens and closes with the elbow movement. With the latch closed, the string is held as the elbow closes. This allows the user to hold a grip while still using the elbow motion. An example where the latch is useful is when the user needs to hold a spoon while eating or a grabber bag while walking. Move the latch slider back to the open position to release the grip.



**Note:** Assembly videos have been prepared and are supplemental or maybe used in place of this document. See at [www.onemonkclapping.com](http://www.onemonkclapping.com)

## Overview

### Phase 1 : Preparation

- Gather Materials
- Identify and Arrange Parts
- Prepare the Printed Parts

### Phase 2 : Hand

- Assembling the Hand
- Stringing the Hand
- Tensioning

### Phase 3: Forearm

- Assembling the Forearm
- Assemble the Cuff

### Phase 4: Putting it all together

- Attaching Hand to Arm
- Assembling the Latch
- Final Tensioning and Knotting

### Phase 5: Final Tests and Fitting

- Straps, Finger Tips, and Palm Grip
- Padding and Fitting

## Gather Materials

### Tools needed for Assembly

We tried to keep the required tools list as short as possible.



Thick **CA glue** such as **Super Glue** for sealing the knots in the fishing line and for assembling the forearm.



**Two part epoxy** is better to assemble arm pieces, if you print the arm in two pieces.



Phillips head **screwdriver** to fit the Sheet metal screws (#4 pan head)



Hexagonal **allen wrench** or Robertson **square head** drive. Depends if you use metal purchased bolts or the 3d printed bolts.

Sharp **knife or scissors** to cut the fishing line.

### Tools Recommended for Assembly

These tools will make your assembly easier and generally result in a better device.



**Small flat file , Small round file , Medium half round file sanding blocks or sandpaper** to clean printed parts.



**Toothpick** or similar disposable fine tip to help spread the glue.



**Utility knife or razor** to help clean parts. A sharp razor is helpful to clean “stringing” off the flexible parts and is very useful to trim the hinges.



A pair of **needle nose pliers** can come in handy, particularly when removing scaffolding.



### **Galvanized Picture Hanging Wire and/or a T handle Ball End Allen Wrench**

Both are very useful in clearing the string holes in the fingers and palm if they close up during 3-D printing. The wire helps pass the string in difficult places. Galvanized wire comes in many sizes at the local hardware store. It is generally stiff enough to force through the holes to open them up.



### **Plasti Dip - clear**

We paint this on the fingertips and the entire palm grip area. This creates a much better grip, then the plastic by itself. It is durable when used with the primer.



A **thread Tap tool** is handy. You can always heat a metal bolt and use it to clean the threads in the plastic. But a tapping tool is made to do this. You likely will only need a few sizes which correspond to the size of your palm and elbow bolts.

## **Parts and Supplies**

These are the materials you will use for each build.



**Dyneema fishing line—80 to 120 test strength** Dyneema is the name of material used in recommended fishing line. For a 100% size hand you will need about 10 feet. Purchase extra to allow for adjustments in tensioning and restringing of the hand. (e.g. brand names include Spiderwire or Hercules).



**Steel rod of 2mm diameter.**

2mm rods are available on Amazon.

<https://www.amazon.com/gp/product/B07D2WSF29>

Typical wire clothes hangers are suitable or hardware stores sell this diameter steel.



We only need a small piece to stabilize the latch mechanism.



**Self Adhesive Padding** in areas intended to be in contact with the recipient's skin (cuff and forearm) is recommended. The thickness of this padding should have already been designated in the OpenSCAD "padding thickness" choice section. If the leather option is chosen for the cuff, no padding is needed in the cuff-but recommended for the socket.

Silicone: <https://amputeestore.com/products/silipos-pressure-relief-gel-padding>

Neoprene:

<https://www.amazon.com/Insulation-Neoprene-Adhesive-Multi-Function-Soundproof/dp/B08R9LF7BX>

or

Neoprene:

<https://www.amazon.com/NATGAI-Neoprene-Adhesive-Multiple-Dimensions/dp/B07VY7SFRP>



**Velcro Straps** —the width of these straps should have been already designated in the OpenSCAD "StrapWidth" choice option. Generally, 8" to 12" long, 1" wide, non-elastic (often sold as cable straps). For small arms you can use 1/2" wide



**Palm bolts Cap Head** It is recommended to use store bought nylon bolts because they are lighter than metal bolts and more durable than the 3-D printed ones.

The bolts need to match the thread size specified in the OpenSCAD. Most hardware stores will carry these bolts.





**Elbow Bolts round head** - A nice smooth round head store bought nylon bolt is the best option. When printing with PETG I have had good luck with the printed elbow bolts. But they can break while tightening which can be frustrating to get cleared out.

The bolts need to match the thread size specified in the OpenSCAD. Most hardware stores will carry these bolts.



**Sheet metal screws (#4 pan head)** of appropriate length for attaching socket wrist to Palm (1 ea.) and for attaching LatchCover to Upper Arm (1 ea.) The #4 size is US standard, outside the US it is a 3mm sheet metal screw.



#### **ThreadLocker Blue Gel**

Once you have assembled everything, it is recommended to put on the bolt threads so they do not work loose. You can use glue if you do not have this.



**Rivets and/or Chicago screw/post** Applicable to the leather option only.

For smaller sizes, we prefer copper rivets. They are small and can be cut to length. The down side is you can crush the plastic if you are not careful.

<https://www.amazon.com/gp/product/B08HMYMFWS>



For larger sizes we prefer **Chicago screws**. They are reusable, but only come in a limited number of lengths and diameters.

Make sure the correct diameter of the shaft is set in the OpenSCAD options



**Leather** We recommend ~4.5 oz leather (~1.8mm thick). It is the more expensive option compared to plastic, but we believe it to be more comfortable than the plastic. You can print the leather template SVG file on paper, and use that to determine how much leather you will need.

A **leather punch** is recommended to make accurate holes in the leather where it will attach to the plastic components. A website of one supplier of cowhide leather is [www.montanaleather.com](http://www.montanaleather.com)

**To determine which size and length to purchase for the elbow and palm bolts, use the 3-D printed bolts and the parameters on OpenSCAD. It is helpful to bring the sample 3-D printed bolts with you to the store. Both metal and nylon bolts can be cut if they are too long and if you have the proper tools.**

## Identify and Arrange Parts



Identify the correct parts for the fingers and palm and arrange them in order.

The finger pieces have numbers embossed in them. These correspond to the finger numbers 1 to 5 starting with the thumb.


The location of the number on the fingers will verify that you have a Right or Left Hand.



## Prepare the Printed Parts

### Thoroughly Remove the Scaffolding

This is an important step. Make sure you remove all scaffolding in the corners and edges.

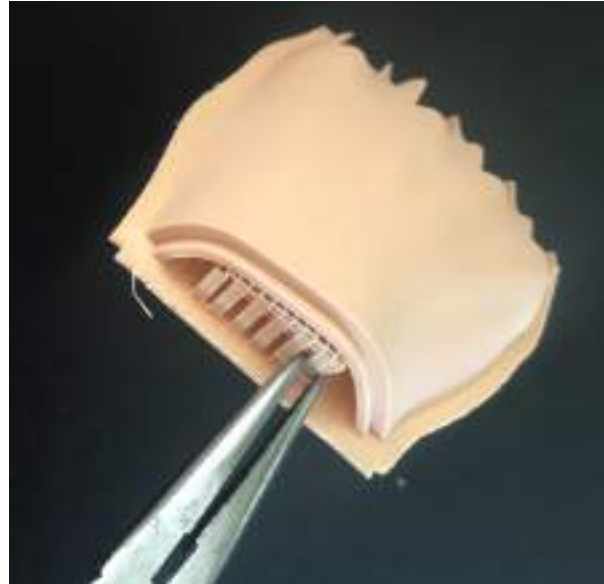

### Sand, File and Smooth Parts

Sand, file and generally finish the parts. If you oriented the finger phalanx as suggested during the printing process, you will need to sand the curved back end of the parts. The added supports and spreading will mean the back is not round.







### Clean string channels

In smaller printed hands the holes for the strings can close during the printing process. Galvanized picture wire is generally stiff enough to force through the holes to open them up.



You can also force the wire through the opening, heat it with a solder iron, and move it around to smooth the inside of the channels.

### Clean the Palm and Elbow threads

Using a metal bolt or a tap tool, try screwing the bolts into just the threaded part. The top of the palmtop and the elbow have threads matching what was set by OpenSCAD.



You may need to heat the bolt to help clean the plastic threads.



### **Clean the wrist bolt threads**

Try screwing the wrist bolt in the forearm. Make sure it is smooth and easy to tighten all the way.

### **Clean wrist button mechanism**

The wrist button is printed with very tight tolerances. You will likely need to sand the ratcheting part of the wrist bolt to a smooth round finish.

Insert the wrist bolt from the inside of the palm. Make sure it rotates easily, yet is snug




Put the Wrist cover over the wrist bolt, **WITHOUT** the WRIST BUTTON. Make sure the wrist bolt can rotate smoothly and easily yet does not rattle. Take the wrist bolt back out after the test.



## B. Hand Assembly

1. Gather all parts of the hand together. Each finger end and phalanx has a number that is imprinted on it, i.e. 1=thumb, 2=index, 3=middle, 4=ring, and 5=pinky. Fit the pieces together—all should fit and move without undo friction that might impede bending. It may be necessary to use a file, utility knife, or sandpaper to facilitate finger movement. The rectangular holes should also be free of burrs or print filament webbing.
2. Insert the TPU hinges in the slots between the finger ends and the finger phalanxes. Insert other TPU hinges in the phalanx top and slide the ends into the slots in the Palm.
3. Stringing the hand is fairly straightforward and much easier than previous Kwawu versions. Figure 1 shows the fingers and the stringing route. Figure 1 along with Figure 2 shows a hand strung appropriately.
4. There is a structure at the end of each finger for tying a knot. Use Super glue to fix the knots from slipping. Note that the long end of the secondary whippetree is toward the thumb and String 1 from the thumb connects to it (Figure 2). String 2 connects index finger end to middle finger end via the thumb end of the primary whippetree. Note the imprinted arrow on the primary whippetree will point toward the thumb. String 3 connects the ring finger end to pinky finger end via the end opposite the thumb end of primary whippetree.
5. Figure 2 shows how the whippetrees are connected to the fingers (Strings 1,2, and 3) and the whippetrees to each other (String 4). String 5 connects to the center hole in the Secondary Whippetree and extends to the Tensioner (Tensioner not shown). String 5 needs to be long enough to reach from the secondary whippetree to the bicep areas of the recipient's arm. It is recommended to attach String 5 later in the assembly process (see Step 7 in Forearm Assembly).

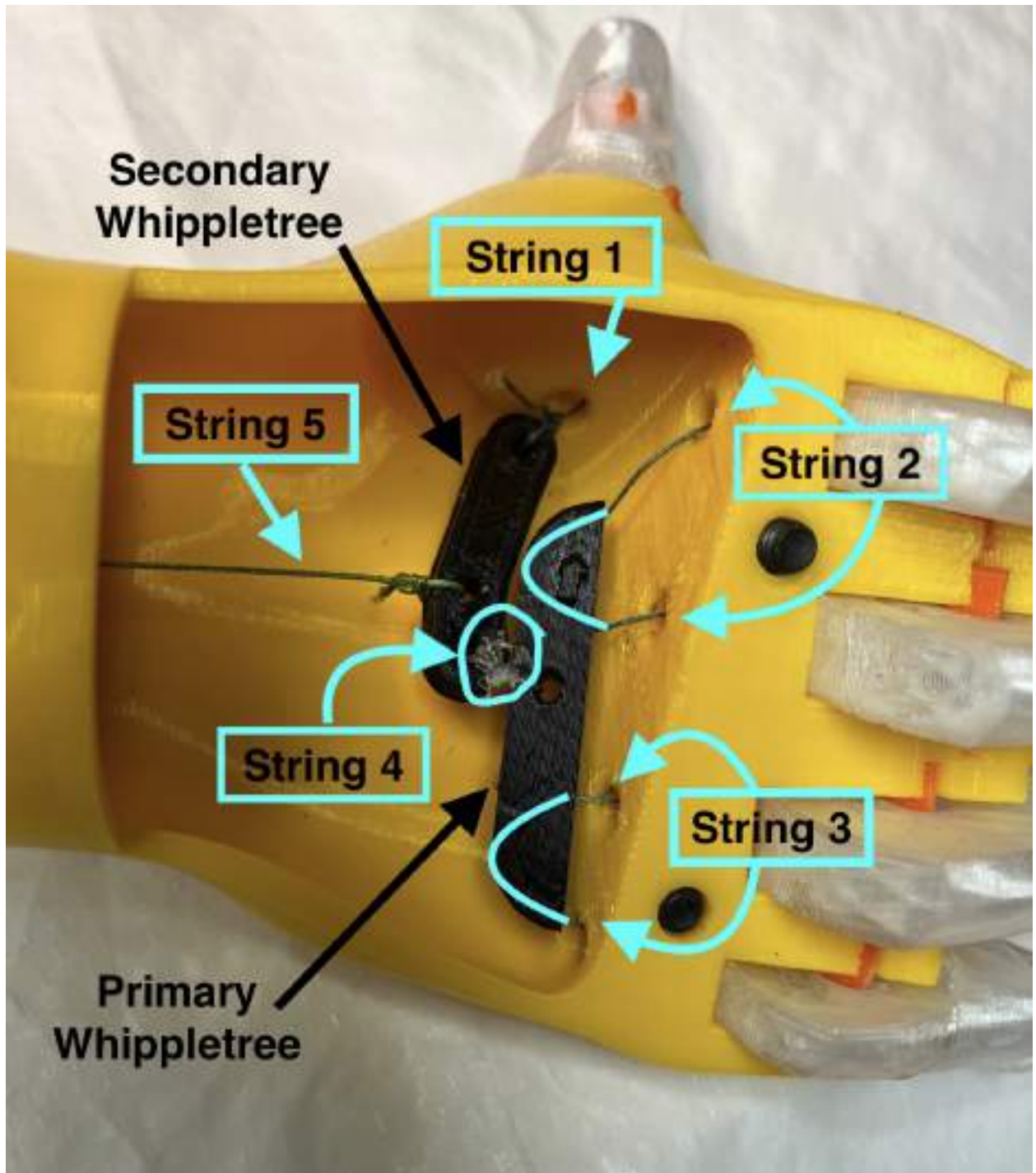

6.

**Figure 2: Stringing the hand: Palm**

7. As shown in Figure 2, with all fingers extended, the the Primary Whippletree should be close to the area in the palm where Strings 2 and 3 enter the palm area. When the Whippletrees are tied together (String 4) and thumb extended, the thumb side of secondary whippletree should be close to the palm area where String 1 enters the palm. This will give the whippletrees adequate room to move when fingers and thumb are clenched.
8. Figure 3 shows components in the wrist. To make sure parts will function properly when assembled, first make sure that the WristButton and WristCover will fit over the WristBolt all the way to the WristBolt pedestal and are able to turn freely. The WristCover was designed to be a tight fit, so it may be necessary to grind/file the ID of the WristCover so that it will fit and turn properly. Also the WristBolt must fit into the hole on the palm top and be able to turn freely. Again, the ID of the Palm top hole may have to be expanded by filing/grinding.




**Figure 3 : Wrist Components**

9. Test the WristBolt to make sure it threads into the ForearmArm. Filing in the threads of the bolt may be required to initially get the WristBolt fully engage in the Forearm threads. The application of a dry lubricant to the threads will help to obtain smooth turning of the bolt into the forearm. After obtaining smooth threading of the bolt into the socket, remove the WristBolt from the forearm.
10. Place the WristButton into the top of Palm cavity. While pushing the button, insert the WristBolt into Palm cavity and push the threads out of the top opening. Then place the WristCover over WristButton such that the wedge shaped structures fit with the recesses in the Palm top.

## C. Forearm Assembly




1. If the forearm was printed in two parts, glue the UpperArm and LowerArm together with two part epoxy.
2. Prepare to assemble the latch mechanism—Latch, LatchCover, Ratchet, and Tensioner onto the UpperArm. See Figure 4.

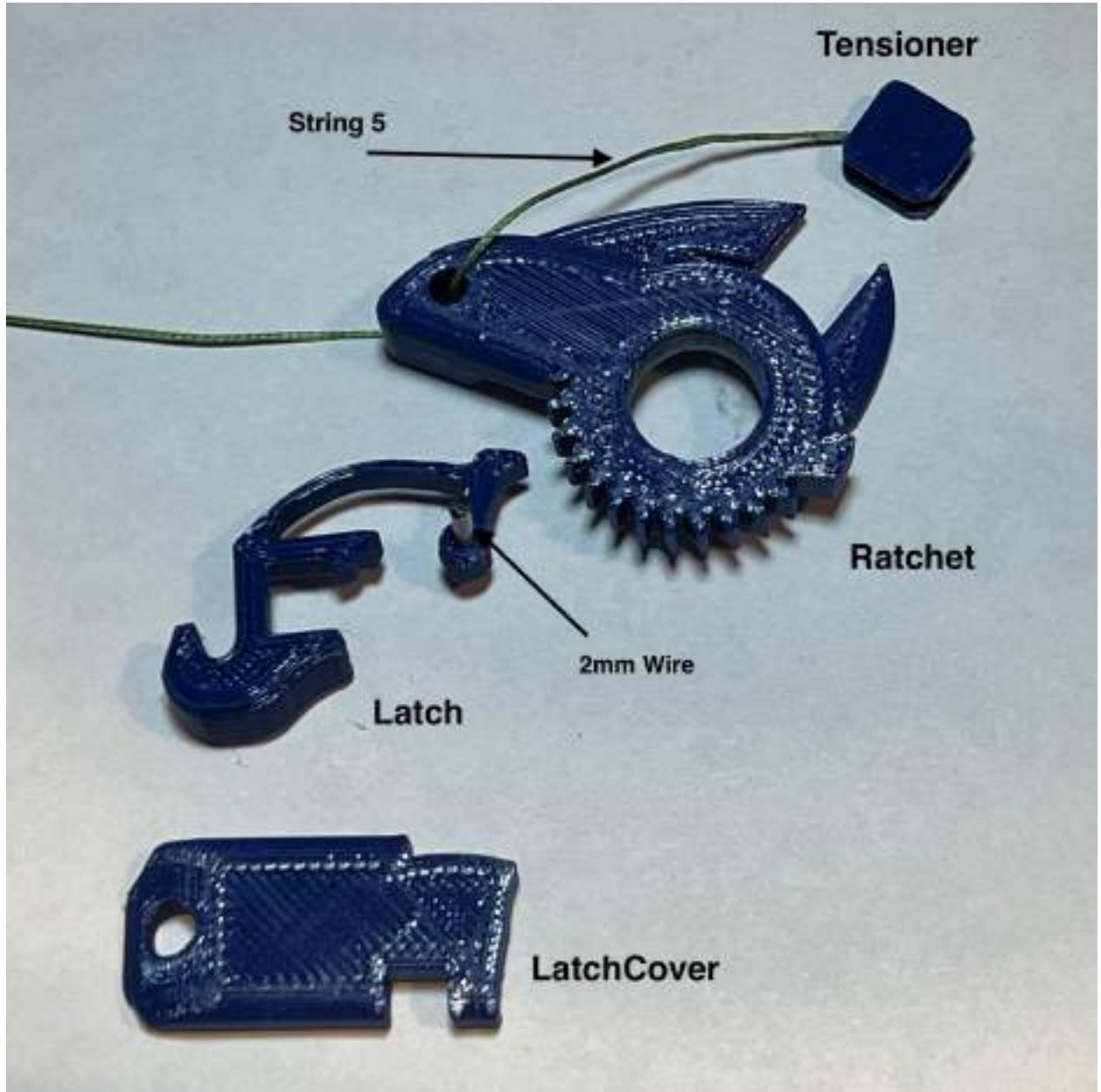

**Figure 4: Latch Mechanism**

3. Cut a piece of 2mm wire/hanger to span from the hole in the UpperArm Latch recess, through the Latch and into the hole in the LatchCover. The hole in the Latch may have to be widened to be able to insert the wire—use appropriate sized drill bit if necessary. Insert the wire in the Latch hole to act as an axel in the Latch and set the Latch with wire into the forearm recess.

4. Place the LatchCover over the Latch making sure the wire is inserted in the hole supplied. Use a #4 pan head sheet metal screw to attach the LatchCover to the Forearm.
5. With the Latch pulled back (off position), place the Ratchet in place. Make sure the Latch will engage in the Ratchet teeth when in the forward (on) position.
6. Insert the Tensioner in the slot in the side of the Ratchet (see Figure 4). It should fit in all four of the possible orientations. Some trimming of the Tensioner may be required in certain orientations.
7. Push one end of String 5 through the Egress Port in the side of the forearm and down through the threaded region. See Figure 5. With the hand still disconnected from the forearm, continue pushing the string through the hole in the WristBolt and into the Palm cavity. Attach the string to the center hole of the Secondary Whippletree as shown in Figure 2.
8. Screw the WristBolt on the Forearm until tight, then back off until you see (through the ForearmWrist Lock hole) the first of the 4 holes that is in the WristBolt. Use a sheet metal screw to secure the Wristbolt to the Forearm. (See Figure 5)

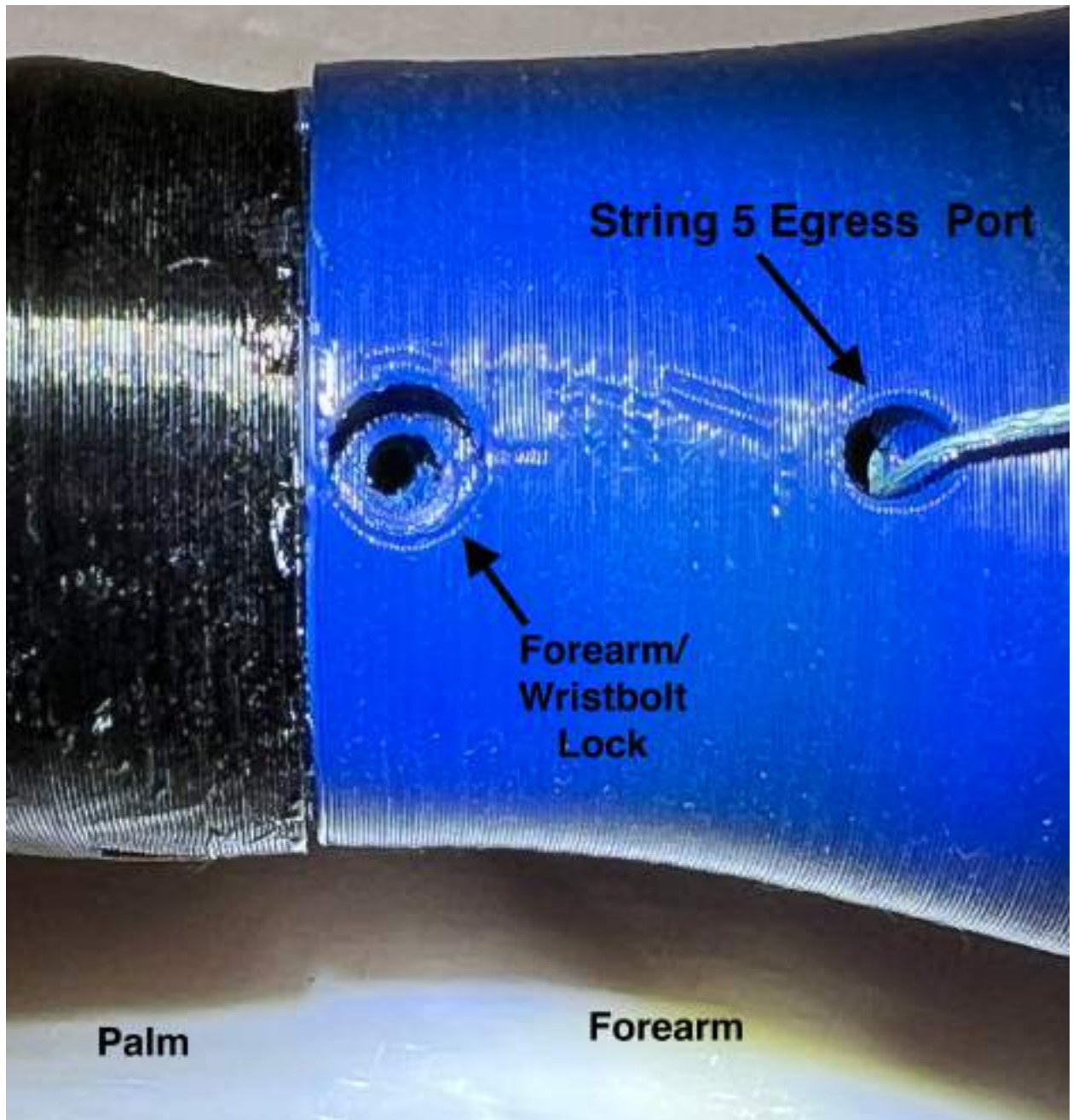

**Figure 5: Forearm to Palm**

9. Cut String 5 such that it extends past the Ratchet on the forearm. Thread the free end of the string up through the hole in the Ratchet as shown in Figure 4, then tie it to the Tensioner through the center hole in its axle. Super glue the both knots in String 5. The PalmCover can now be attached over the palm cavity by fastening it with 2 bolts selected earlier.



## D. Cuff Assembly

1. Plastic Cuff: The Plastic Cuff is ready to connect with the Forearm as printed
2. Leather Cuff:
  - a. Lay the leather (cut previously) flat with the top smooth surface that will eventually be in the contact with the skin of the recipient. Place the printed Cuff1 under the left side of the leather and Cuff2 under the right side. Verify the holes in the cuffs line up with the holes in the leather.
  - b. Attach the leather to the Cuff1 and Cuff2 printed parts using Rivets or Chicago type fasteners. Push the fasteners through the leather and plastic and seal from the plastic side.
3. Both Plastic and Leather Cuffs
  - a. Using bolts previously selected, attach Cuff1 to the Ratchet/Latch side of the forearm and attach Cuff2 to the opposite side. It may be helpful to use appropriate size taps to clean out the printed threads.
  - b. Flex the Cuff and verify finger movement. If String 5 is not tight enough, remove the Tensioner from the Ratchet slot and wind the string tighter and place it back in the Ratchet slot.
4. Velcro Strap Install
  - a. Use a strap having the width as chosen initially and thread it through the gaps between the plastic or leather and through the attached loop at the velcro end. The strap will allow the cuff to be tightened as necessary

## E. Padding Install

1. Plastic Cuff:
  - a. Silicone Option: Cut the silicone in a pattern to fit inside the plastic grid on the cuff. Glue the silicone to the grid with a flexible adhesive.
  - b. Neoprene Option: Cut the self adhesive neoprene to fit inside the plastic grid on the Cuff. Peel the backing and stick it onto the grid.
2. Leather Cuff—the leather Cuff does not need padding
3. Socket
  - a. The padding needs to be just long enough to cover the recipient's residual but shorter as to not cover any existing vent holes. Make a paper template of the inside of the socket and use it as a pattern for silicone or neoprene padding.
  - b. Silicone Option: Cut the silicone to the pattern and dry fit it inside. the socket-some trimming may be required. Glue the silicone to the inside to the socket with a suitable adhesive.
  - c. Neoprene Option: Cut the neoprene to the pattern and dry fit it inside the socket-some trimming may be required. Remove the backing paper and self adhere the neoprene inside the socket.

## F. Final Comments

1. Plasti Dip is commonly used on fingers and palms to enhance the ability of the Kwawu 3 to pick and hold things because it adds tackiness to the plastic fingers. Plasti Dip Primer is recommended prior to painting with the Plasti Dip. The joints and strings in the fingers should not be exposed to Plasti Dip.
2. Lee Tippy nail's Gel Finger Tips can be purchased and placed over the ends of the plastic fingers to enhance gripability.

## Final Testing

1. Exercise the wrist lock mechanism to make sure the hand can be turned 360 degrees and can be locked at any selected location.
2. Exercise the latch mechanism. In the engaged position, the latch should lock the fingers in position depending on the extent of elbow rotation.
3. With the latch in the engaged position, fingers should maintain their clutched position even when the elbow is straightened,